


Instructions for use



CE
0297

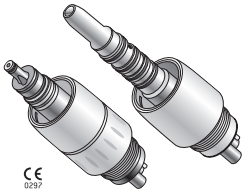

rote quick

**Quick couplings RQ-03, RQ-04,
RQ-14, RQ-24, RQ-34
RA-24, RA-25**

Contents

W&H Symbols	4
1. Introduction	6
2. Safety notes	10
3. Product description	14
4. Operation	16
Connecting the Quick coupling	16
Disconnecting the Quick coupling	17
Changing air driven products	19
Test run	22
5. Hygiene and maintenance	23
General notes	23
Limitations on processing	25

Initial treatment at the point of use.....	26
Manual cleaning.....	27
Manual disinfection	28
Inspection, Maintenance and Testing	29
Packaging.....	30
Sterilization	31
Storage	34
6. Maintenance	35
7. Servicing.....	40
8. W&H Accessories and spare parts	41
9. Technical data	42
10. Disposal	44
Explanation of warranty terms.....	45
Authorized W&H service partner	47

W&H Symbols

in the Instructions for use



WARNING!
(Risk of injury)



ATTENTION!
(to prevent damage
occurring)



General explanations,
without risk to
persons or objects



Do not dispose of
with domestic waste

W&H Symbols

on the medical device/packaging




Medical Device





Type B applied part
(not suitable for intracardiac application)

W&H Symbols


on the medical device/packaging


 CE marking
with identification number
of the Notified Body


 DataMatrix Code
for product information including
UDI (Unique Device Identification)


 HIBC Data structure in accordance
with Health Industry Bar Code


 REF Catalogue number

 UL Component Recognition Mark
indicates compliance with Canadian
and U.S. requirements

 135°C
Sterilizable up to the
stated temperature

 SN Serial number

 **R_x**_{only} Caution! Federal law restricts this device to sale by or on the order of a
dentist, physician, veterinarian or with the descriptive designation of any
other practitioner licensed by the law of the State in which the practitioner
practices to use or order the use of the device.

 Date of manufacture

1. Introduction

Customer satisfaction is the main priority under the W&H quality policy. This medical device has been developed, manufactured and subjected to final inspection according to legal regulations, quality norms and industry standards.

For your safety and for the safety of your patients

Prior to initial use please read the Instructions for use.

These explain how to use your medical device and guarantee a smooth and efficient operation.



Observe the safety notes.

Intended use

Connector for media transfer (air, water, electricity and light) between the supply hose and the dental unit and air driven motors



Misuse may damage the medical device and hence cause risks and hazards for patient, user and third parties.

Qualifications of the user

We have based our development and design of the medical device on the dentists, dental hygienists, dental employees (prophylaxis) and dental assistants target group



0297

Production according to EU Directive

The medical device meets the requirements of Directive 93/42/EEC.

Responsibility of the manufacturer

The manufacturer can only accept responsibility for the safety, reliability and performance of the medical device when it is used in compliance with the following directions:

- > The medical device must be used in accordance with these Instructions for use.
- > Only the components approved by the manufacturer may be replaced (e.g. o-rings).
- > Modifications or repairs must only be undertaken by an authorized W&H service partner. (see page 47)
- > Correct the malfunction as described in the Instructions for use.

Skilled application

The medical device is intended only for skilled application according to the intended use as well as in compliance with the valid health and safety at work regulations, the valid accident prevention regulations and in compliance with these Instructions for use.

The medical device should be prepared for use and maintained by staff who have been trained in procedures for infection control, personal safety and patient safety.

Improper use, (e.g., through poor hygiene and maintenance), non-compliance with our instructions or the use of accessories and spare parts which are not approved by W&H, invalidates all claims under warranty and any other claims.



Any serious incident that has occurred in relation to the medical device should be reported to the manufacturer and the competent authority!

2. Safety notes



- > Before using the medical device for the first time, store it at room temperature for 24 hours.
- > The operation of the medical device is permitted only on dental units which correspond to the standards IEC 60601-1 (EN 60601-1) and IEC 60601-1-2 (EN 60601-1-2).

The power supply unit for the dental unit must satisfy the following requirements to be guaranteed by the system assembler (Relates to externally electrically supplied couplings and applied parts):

- > Double insulation for the highest expected supply voltage must be provided between the primary and secondary power circuits.
- > Double insulation for the highest expected secondary voltage must be provided between the secondary voltage and protective earth (PE).



- > The secondary circuits must be galvanically isolated from each other.
 - > The secondary circuits must be protected against short-circuiting and overloading. (maximum 700 mA only for RA-24)
 - > The leakage currents of and between the applied parts must be kept.
 - > The secondary voltage to supply this medical device must be limited to a maximum of 4.2 V AC or 6 V DC. (3.3 V AC/DC only for RA-24)
-
- > Never touch the patient and the electrical contacts on the medical device simultaneously.
 - > Use only the supply hoses as specified by EN ISO 9168.
 - > Always ensure the correct operating conditions and cooling function.
 - > Check the medical device for damage and loose parts before using.
 - > Do not operate the medical device if it is damaged.
 - > Perform a test run each time before using.
 - > Do not look directly into the light source.
 - > Do not use the medical device as a light probe.



- > The Quick coupling is a functional part of the supply hose and should therefore also be seen as an extension to it during reprocessing. It is imperative to comply with the concentrations and exposure times specified by the manufacturer of the treatment water decontamination system, as well as its handling. (except RQ-14, RQ-34)

If the Quick coupling is processed separately from the supply hose, you can refer to the information in the chapter “Hygiene and maintenance” as per ISO 17664 from the manufacturer of the quick coupling.



The medical device is not approved for operation in potentially explosive atmospheres (only for RA-24).

Hygiene and maintenance prior to initial use




- > The medical device is sealed in PE film and not sterilized when delivered.
- > The PE film and the packaging are non-sterilizable.



- > Clean and disinfect the medical device.
- > Sterilize the medical device (only for sterilizable types).

3. Product description

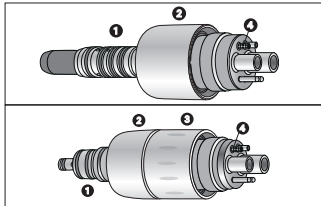

- ❶ O-rings
- ❷ Retention sleeve
- ❸ Spray regulation ring (for RQ-14, RQ-34)
- ❹ Water filter with resuction stop



All Quick couplings are equipped with a non-retraction valve which prevents contaminated cooling water from being sucked back into the turbine and the supply hose.

The non-retraction valve is integrated in the cooling water supply system.

In the event of blocked or incorrectly routed cooling water lines, please contact an authorized W&H service partner (see page 47).




The cooling water lines must not be cleaned with sharp objects!

(This could damage the sealing element and prevent the non-retraction valve from working.)

4. Operation

Connecting the Quick coupling

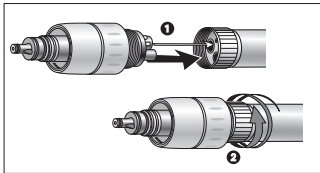

- 1 Connect the Quick coupling to the supply hose.
- 2 Firmly tighten the union nut of the supply hose by hand in a clockwise direction to ensure there are no leaks.




Verify full engagement.
Check the leak tightness.



This method of assembly provides a connection for drive air, chip air, return air, water.
Possible variants with electricity: RQ-24, RQ-34 or with light: RA-24, RA-25.

Operation

Disconnecting the Quick coupling

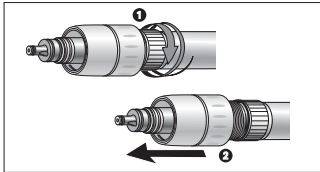

- 1 Unscrew the union nut of the supply hose by hand in an anticlockwise direction.
- 2 Carefully remove the Quick coupling from the supply hose.


Checking the Quick coupling for leaks



- > Push an appropriate air driven medical product onto the Quick coupling.
- > Activate the medical device, or if possible just activate the spray water only.
- > No water should leak between the Quick coupling and the air driven medical product, and the coupling and the supply hose.



- > Always follow recommendations made by the manufacturer of air driven products.
- > Only connect air driven product with appropriate connection to the Quick coupling.
- > The user accepts sole responsibility if other air driven products are used. We accept no liability in such cases.

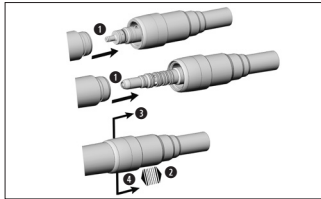

Do not assemble or remove the medical device during the operation.

- 1 Attach the air driven instrument onto the Quick coupling until it audibly engages.



- 2 Verify full engagement

- 3 Pull the retention sleeve of the Quick coupling back.
- 4 To disconnect the air driven instrument, pull it off in axial direction.

Regulating the spray water

The water quantity is regulated directly via the spray control ring of the medical device (RQ-14, RQ-34) or via the adjustable valves on the unit (RQ-03, RQ-04, RQ-24, RA-24, RA-25).

Minimum spray flow: Line marking on the spray adjuster ring and point on the retention sleeve correspond with each other.

Maximum spray flow: Turn the spray adjuster ring to the right or left.



The water flow can only be varied using this spray adjuster ring. It cannot be stopped completely.

Test run








Do not hold the medical device at eye level.

- > Connect the Quick coupling to the supply hose.
- > Attach the air driven instrument onto the Quick coupling until it audibly engages.
- > Start the air driven product.



In the event of operating malfunctions (e.g. vibrations, unusual noise, overheating, coolant failure or leakage) **stop the medical device immediately** and contact an authorized W&H service partner.

-  > Follow your local and national laws, directives, standards and guidelines for cleaning, disinfection and sterilization.
-  > The information on the validated reprocessing procedures serves as an example of an ISO 17664 compliant processing of the medical device.
-  > Wear protective clothing, safety glasses, face mask and gloves.
-  > Couplings are considered an extension of the tubing. Clean and disinfect without disconnecting from the tubing using an intermediate, hospital grade disinfectant after each patient.
-  > Use only oil-free, filtered compressed air with a maximum operating pressure of 3 bar (43.5 psi) for manual drying.


Cleaning agents and disinfectants




- > Read the notes, follow the instructions and heed the warnings provided by the manufacturers of cleaning agents and/or disinfectants.
- > Use only detergents which are intended for cleaning and/or disinfecting medical devices made of metal and plastic.
- > It is imperative to comply with the concentrations and exposure times specified by the manufacturer of the disinfectant.
- > Use disinfectants which have been tested and found effective by the Verbund für Angewandte Hygiene e. V. (VAH = Association for Applied Hygiene), the Österreichische Gesellschaft für Hygiene, Mikrobiologie und Präventivmedizin (ÖGHMP = Austrian Society for Hygiene, Microbiology and Preventive Medicine), the Food and Drug Administration (FDA) and the U.S. Environmental Protection Agency (EPA).





The user is responsible for validating its process if the specified cleaning agents and disinfectants are not available.


- 
- The product lifetime and the medical device's ability to operate correctly are mainly determined by mechanical stress during use and chemical influences due to processing.
- > Send worn or damaged medical devices and/or medical devices with material changes to an authorized W&H service partner.


Processing cycles

- 
- > In the case of wipe disinfection, the use of the medical device is guaranteed without restriction until a functional or material limitation is recognizable.
 - > The use of the medical device is guaranteed until at least 500 reprocessing cycles.

-  Clean the medical device immediately after every treatment.
 - > Wipe the entire surface of the instrument with disinfectant.

-  > If the Quick coupling remains on the supply hose, follow the instructions of the unit manufacturer.

-  > If the Quick coupling is to be reprocessed, remove it from the supply hose and follow the steps on page 27 - 34.

-  Note that the disinfectant used during pre-treatment is only for personal protection and cannot replace the disinfectant step after cleaning.



Do not place the medical device in liquid disinfectant or in an ultrasonic bath.

- > Clean the medical device under running tap water (< 35°C / < 95°F).
- > Rinse and brush off all internal and external surfaces.
- > Remove liquid residues using compressed air.



> W&H recommends wipe-down disinfection.



Evidence of the medical device's basic suitability for effective manual disinfection was provided by an independent test laboratory using the disinfectants "mikrozid® AF wipes" (Schülke & Mayr GmbH, Norderstedt) and "CaviWipes™" (Metrex).



The medical device is not approved for automated cleaning and disinfection.

Inspection



- > Check the medical device after cleaning and disinfection for damage, visible residual soiling and surface changes.
- > Reprocess any medical devices that are still soiled.
- > Sterilize the medical device following cleaning and disinfection.



Pack the medical device in sterilization packages that meet the following requirements:

- > The sterilization package must meet the applicable standards in respect of quality and use and must be suitable for the sterilization method.
- > The sterilization package must be large enough for the sterilization goods.
- > The filled sterilization package must not be under tension.



W&H recommends sterilization according to EN 13060, EN 285 or ANSI/AAMI ST55.



- > Read the notes, follow the instructions and heed the warnings provided by the manufacturers of steam sterilizers.
- > The program selected must be suitable for the medical device.

Recommended sterilization procedures

- > “Dynamic-air-removal prevacuum cycle” (type B) / “Steam-flush pressure-pulse cycle” (type S)*/**
134°C (273°F) for at least 3 minutes, 132°C (270°F) for at least 4 minutes
- > “Gravity-displacement cycle” (type N)**
121°C (250°F) for at least 30 minutes
Maximum sterilization temperature 135°C (275°F)



Evidence of the medical device's basic suitability for effective sterilization was provided by an independent test laboratory using the LISA 517 B17L* steam sterilizer (W&H Sterilization S.r.l., Brusaporto (BG)), the Systec VE-150* steam sterilizer (Systec) and the CertoClav MultiControl MC2-S09S273** steam sterilizer (CertoClav GmbH, Traun).

- > "Dynamic-air-removal prevacuum cycle" (type B): 134°C (273°F) – 3 minutes*, 132°C (270°F) – 4 minutes*/**
- > "Steam-flush pressure-pulse cycle" (type S): 134°C (273°F) – 3 minutes*, 132°C (270°F) – 4 minutes*/**
- > "Gravity-displacement cycle" (type N): 121°C (250°F) – 30 minutes**

Drying time:

- > "Dynamic-air-removal prevacuum cycle" (type B): 132°C (270°F) – 30 minutes**
- > "Steam-flush pressure-pulse cycle" (type S): 132°C (270°F) – 30 minutes**
- > "Gravity-displacement cycle" (type N): 121°C (250°F) – 30 minutes**

* EN 13060, EN 285, ISO 17665


** ANSI/AAMI ST55 , ANSI/AAMI ST79



- > Store sterile goods dust-free and dry.
- > The shelf life of the sterile goods depends on the storage conditions and type of packaging.




The Quick coupling may be stored on the supply hose.

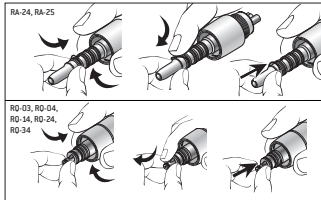

- > Replace damaged or leaking o-rings immediately.
- > Always replace all o-rings.
- > Do not use sharp tools.

- 1 Press the o-ring firmly together with your thumb and index finger until it forms a loop.
- 2 Pull off the o-ring.
- 3 Slide the new o-ring back on.

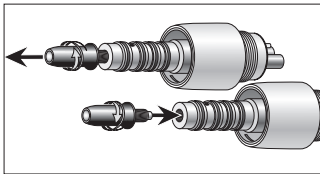


Replacing the bulb (RA-24)

- > Screw off the bulb counterclockwise by hand and replace it with a new one.




Verify full engagement.

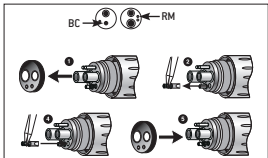

Changing the water filter (BC/RM)

- 1 Remove the seal.
- 2 Pull out the water filter with tweezers.
- 3 Clean the water filter (see page 38).
- 4 Carefully insert the water filter.
- 5 Slide on the seal.

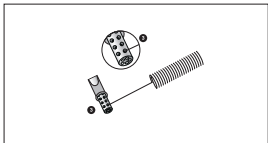

Clean the water filter (BC/RM)

- 3 Use the nozzle cleaner carefully to remove dirt and deposits from the water filter outlets.




The water filter can be cleaned in an ultrasonic bath.



If it proves impossible to correct the malfunction, please contact an authorized W&H service partner.



> Repeat the complete hygiene and maintenance process.

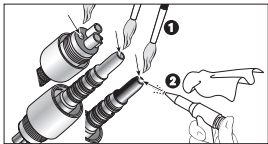

Cleaning of the light source (RA-24, RA-25)



Avoid scratching the light source!

- ❶ Wash the light source with cleaning fluid and a soft cloth.
- ❷ Blow the light source dry using compressed air or dry it carefully with a soft cloth.



- > Carry out a visual inspection after each cleaning process.
- > Do not use the medical device if the optic outlet is damaged and contact an authorized W&H service partner.

In the event of operating malfunctions immediately contact an authorized W&H service partner.
Repairs must only be undertaken by an authorized W&H service partner.



> Ensure that the medical device has been completely processed before returning it.

8. W&H Accessories and spare parts



Use only original W&H accessories and spare parts or accessories approved by W&H.

Suppliers: W&H partners

02060100	O-ring, large
02060200	O-ring, small
02207300	RM seal
01000700	BC seal (RQ-03)
04697000	Back-pressure valve
04697100	Water tube
03523000	Spare bulb (RA-24)

9. Technical data

Quick coupling	RQ-03	RQ-04	RQ-14	RQ-24	RQ-34	RA-24	RA-25
Connection according to standard EN ISO 9168:2009	Type 1: "Borden"	Type 2: "Ritter Midwest (4 hole)"		Type 3: "Ritter Midwest (6 hole)"			Type 2: "Ritter Midwest (US)"
Connection medical device	"Roto Quick" products					"LW" products	
Water quantity adjustable			x		x		
Light						x Halogen lamp	x Glass rod
Electrical contacts (for power transmission to the medical device)				x	x		
Recommended voltage range V DC or V AC $\pm 0,1$				3.2	3.2	3.2	
Sterilizable	x	x	x	x	x	x	x



Temperature information

Temperature of the medical device on the operator side: maximum 55°C (131°F)

Ambient conditions

Temperature during storage and transport:	-40 °C to +70 °C (-40 °F to +158 °F)
Humidity during storage and transport:	8 % to 80 % (relative), non-condensing
Temperature during operation:	+10 °C to +35 °C (+50 °F to +86 °F)
Humidity during operation:	15 % to 80 % (relative), non-condensing

10. Disposal



Ensure that the parts are not contaminated on disposal.



Follow your local and national laws, directives, standards and guidelines for disposal

> Medical device

> Packaging

Letter of indemnity

This medical device has been manufactured with great care by highly qualified specialists. A wide variety of tests and controls guarantee faultless operation. Please note that claims under warranty can only be validated when all the directions in the instructions for use have been followed.

As manufacturer, W&H is liable for material or manufacturing defects within a warranty period of 24 months from the date of purchase.

We accept no responsibility for damage caused by incorrect handling or by repairs carried out by third parties not authorized to do so by W&H!

Claims under warranty – accompanied by proof of purchase – must be sent to the vendor or to an authorized W&H service partner. The provision of service under warranty extends neither the warranty period nor any other guarantee period.

24 months warranty

Authorized W&H service partner

Find your nearest authorized W&H service partner at <http://wh.com>
Simply go to the menu option "Service" for full details.

Or simply scan the QR code.



Manufacturer

W&H Dentalwerk Bürmoos GmbH
Ignaz-Glaser-Straße 53, 5111 Bürmoos, Austria

t + 43 6274 6236-0, f + 43 6274 6236-55
office@wh.com wh.com

Form-Nr. 50610 AEN
Rev. 009 / 10.02.2021
Subject to alterations