

Technical data

CE
0297

Type:	PB-200 LED	PB-200
Handpiece:	5x LED ring	without light
Supply hose with coupling:	1.8 m	1.8 m
Handpiece weight:	52.5 g	52.5 g
Coolant hose:	2 m	2 m
Coolant volume:	approx. 50 ml/min	approx. 50 ml/min
Mains voltage:	100 – 240 V, 50/60 Hz	100 – 240 V, 50/60 Hz

Equipment supplied Pyon 2 / PB-200 LED:

- › Control unit (100 – 230 V), including 1,8 m supply hose with coupling
- › Handpiece PB-3 LED
- › Tip 1U, 3U, 1P with tip changer and tip measuring card
- › Foot control C-NF
- › Switch-mode power supply
- › Coolant hose and coolant filter
- › O-Ring for hose coupling
- › Nozzle cleaner

Equipment supplied Pyon 2 / PB-200:

- › Control unit (100 – 230 V), including 1,8 m supply hose with coupling
- › Handpiece PB-3
- › Tip 1U, 3U with tip changer and tip measuring card
- › Foot control C-NF
- › Switch-mode power supply
- › Coolant hose and coolant filter
- › O-Ring for hose coupling
- › Nozzle cleaner

W&H Dentalwerk Bürmoos GmbH
Ignaz-Glaser-Straße 53, Postfach 1
5111 Bürmoos, **Austria**
t +43 6274 6236-0
f +43 6274 6236-55
office@wh.com
wh.com

#whdentalwerk
f @ in v
video.wh.com

Now at
your dealer
or **wh.com**



You won't believe your eyes



Pictures are for illustrative purposes only. Additional equipment and accessories shown are not included as standard.

20311 AEN Rev. 006 / 05.09.2019
Subject to alterations

pyon 2 led

Easy to use, top results

Pyon 2: the LED piezo scaler for scaling, periodontology, endodontics and micropreparation. With Pyon 2 you can be sure of achieving optimal results in all applications, as the LED ring in the handpiece allows thorough, precise working. This is helpful both for removing interdental calculus and expanding the root canal entrance.

5 small LEDs

One enormous difference in the application.



Large tip range

Best results in endodontics, prophylaxis, periodontology and micropreparation.



Reliability

Pyon 2 complies with the requirements of the IEC 60601-1 3rd edition.



PB-200 LED

Five LEDs: Robust and sterilizable light sources guarantee optimal illumination. Significantly longer life than conventional light sources.



Optimal hygiene

Thermo washer disinfectable and sterilizable: Both handpieces can be easily dismantled for better cleaning.



Easy to operate

thanks to the stepless power regulation.




Ergonomic design

The shape of the handpiece and its light weight allow you to work for longer without tiring.



Piezo built-in kit

For installation in your existing dental unit.

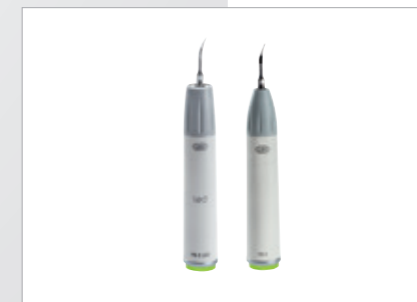

Two versions

The Pyon 2 handpiece is available with a 5x LED ring or without light.